
PROTAC®

INTEGRATING SENSES

Protac® MyRest

- DA Brugsvejledning
- EN Operations manual
- DE Gebrauchsanweisung
- FI Käyttöohje
- FR Mode d'emploi
- NL Gebruiksaanwijzing
- SV Bruksanvisning
- ES Instrucciones de uso
- IT Istruzioni

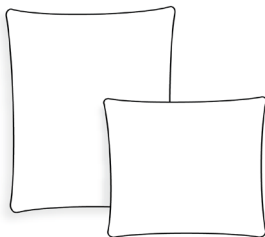Dansk	5
English.....	10
Deutsch	15
Suomi	20
Français	25
Nederlands.....	30
Svenska	35
Español	40
Italiano.....	45

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English

Marking	11
Incidents - adverse events	11
Traceability	11
Key to symbols - product label	11
Safety - general precautions	11
Use of the product	12
Use environment	13
Indications	13
Contra-indications	13
Maintenance	14
Cleaning	14
Key to symbols - care label	14
Service life and warranty	14
Disposal	14

EN

You have chosen a Protac quality product. Protac guarantees that this product is of original Danish design and developed in consultation with experts.

Since 1994, Protac has developed, manufactured and sold sensory-stimulating aids to the healthcare sector in Denmark. Today, Protac products are sold worldwide, especially in Europe, Asia, Australia, New Zealand, Canada and the US.

The Protac products have been developed on the basis of the theories of sensory integration, and they are supported by extensive research, documentation and case studies.

Sensory integration forms the basis in the development of all Protac products. The theory of sensory integration originated in the 1970s in the US, and was developed by A. Jean Ayres, an occupational therapist and psychologist. The theory is based on how the brain processes sensory impression and how the various sensory stimuli can be influenced by applying pressure to the tactile sense (the sense of touch) and the proprioceptive sense (the sense of movement) in muscles and joints. The pressure and movement of the balls inside the product affect these senses in particular.

There is a strong focus on safety and environmental protection during the manufacturing process.

Further information can be found at www.protac.dk.

Marking

The products are CE-marked in accordance with Regulation (EU) 2017/745 on medical devices.

EU declaration of conformity can be found at www.protac.dk

Incidents - adverse events

- Incidents related to the use of Protac® MyRest must be reported to Protac.
- Adverse events related to the use of Protac® MyRest must be reported to the National Competent Authority and to Protac.

Traceability

Production month, year and batch number can be found on the product label.

Do not remove the product label from the product.

EN

Key to symbols - product label

	Instructions for use		Maximum user weight
	Medical device		Article number
	Lot number		For indoor use
	Manufacturer		Production month/year
	CE marked in accordance with Regulation (EU) 2017/745 on medical devices		
	Warning: See under Safety		


Safety - general precautions

	Not for unsupervised use by users unable to care for themselves
	Avoid smoking next to the product. Dropped embers or cigarettes can cause burns.
	Never leave switched on electrical appliances that emit heat on or near the product.
	Do not use near open fire.
	Do not use other balls in the product than those supplied by the manufacturer.


Use of the product

Protac® MyRest is supplied ready for use. To be used with cover (included).


The body pillow can be used to alleviate restless legs, while sleeping/resting on one side.

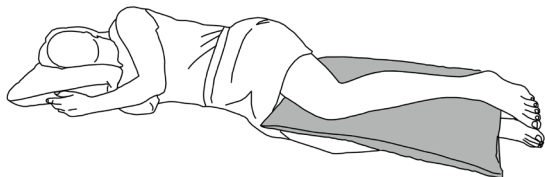The body pillow can also be used while lying on the back.

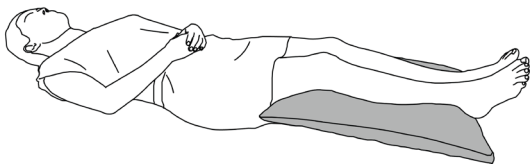Two body pillows can be used at the same time to create a comfortable seating or lying position.

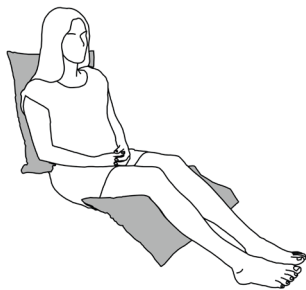Use environment

Protac® MyRest is for indoor use. The body pillow can be used both in institutions and private homes.

Indications

Protac® MyRest is a sensory-stimulating aid intended to alleviate sleeping difficulties caused by restless legs and/or growing pains. Furthermore, it can be used for persons (children, adults, elderly) with sensory impairment.

Protac® MyRest contains balls which move and apply dynamic pressure on to the body, thereby stimulating the muscle-joint sense and the tactile sense. This increases body awareness and calms the brain and body, and the product thereby alleviates physical and psychological unrest.

Contra-indications

Persons with pressure ulcers or persons at risk of developing pressure ulcers are not to use the product.

Using Protac® MyRest as a head pillow might in some cases cause neck soreness.

Maintenance

It is recommended to fluff the pillow on a regular basis.

Cleaning

Protac® MyRest can be washed in a regular household washing machine. The frequency of washing depends on use. When washing, please follow the instructions on the care label. Washing instructions can also be found on pages 50-51 of this manual.

Protac does not recommend the use of stain remover.

Key to symbols - care label

	Maximum washing temperature		Do not iron
	Can be tumble dried		Do not dryclean
	Do not tumble dry		Do not bleach

Service life and warranty

The expected service life of the product is 5 years. This is based on correct use, cleaning, washing and maintenance.

Products used beyond the expected service life must regularly undergo a maintenance inspection which implies checking for wear and tear in fabric and seams.

The product comes with a 2-year warranty, see page 2.

The warranty will be invalid in case of:

- overload or excessive use of the product
- lack of or incorrect maintenance of the product

Disposal

The product can be broken down into the following components for disposal:

- The plastic balls can be disposed of as hard plastic.
- The fabric can be disposed of as textiles.

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Protac A/S