
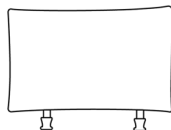
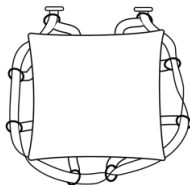
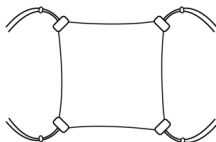
PROTAC®

INTEGRATING SENSES

Protac Kuglepuden® Protac Ball Cushion®

- DA Brugsvejledning
- EN Operations manual
- DE Gebrauchsanweisung
- FI Käyttöohje
- FR Mode d'emploi
- NL Gebruiksaanwijzing
- SV Bruksanvisning
- ES Instrucciones de uso
- IT Istruzioni

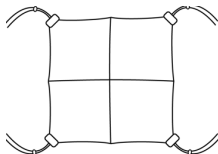Dansk	6
English.....	11
Deutsch	16
Suomi	21
Français	26
Nederlands.....	31
Svenska	36
Español	41
Italiano.....	46

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English

Marking	12
Incidents - adverse events	12
Traceability	12
Key to symbols - product label	12
Safety - general precautions	12
Use of the product	13
Use environment	14
Indications	14
Contra-indications	14
Maintenance	15
Cleaning	15
Key to symbols - care label	15
Service life and warranty	15
Disposal	15

EN

You have chosen a Protac quality product. Protac guarantees that this product is of original Danish design and developed in consultation with experts.

Since 1994, Protac has developed, manufactured and sold sensory-stimulating aids to the healthcare sector in Denmark.

Today, Protac products are sold worldwide.

All Protac products have been developed on the basis of the theories of sensory integration, and they are supported by extensive research, documentation and case studies.

The theory of sensory integration by A. Jean Ayres forms the basis of the development of all Protac products. The theory is based on how the brain processes sensory impression and how stimulation of the muscle-joint-sense and the deep tactile sense has a calming and organising effect on the sensory system.

All Protac products contain balls, and it is the movement of the balls that provide the dynamic deep-point pressure on the body and stimulate the receptors in the deep tactile sense and the muscle-joint sense. This enhances the body awareness, alleviates sensory disorder and has a calming effect on both body and mind.

Further information can be found at www.protac.dk.

Marking

The products are CE-marked in accordance with Regulation (EU) 2017/745 on medical devices.

EU declaration of conformity can be found at www.protac.dk

Incidents - adverse events

- Incidents related to the use of Protac Ball Cushion® must be reported to Protac.
- Adverse events related to the use of Protac Ball Cushion® must be reported to the National Competent Authority and to Protac.

Traceability

EN Production month, year and batch number can be found on the product label.
Do not remove the product label from the product.

Key to symbols - product label

	Instructions for use		Maximum user weight
	Medical device		Article number
	Lot number		For indoor use
	Manufacturer		Production month/year
	CE marked in accordance with Regulation (EU) 2017/745 on medical devices		
	Warning: See under Safety		

Safety - general precautions

	The cushion must be securely fastened to the chair.
	Do not use cushions with cord locks for children under the age of 3 years.
	Avoid smoking next to the product. Dropped embers or cigarettes can cause burns.
	Never leave switched on electrical appliances that emit heat on or near the product.
	Do not use near open fire.
	Do not use other balls in the product than those supplied by the manufacturer.

Use of the product

Protac Ball Cushion® is supplied ready for use.

It is important that the cushion is placed on a rigid seat surface to gain full benefit from the sensory stimulation.



There are different attachment methods for different types of chairs.



On a chair with a single column, the cushion is attached by pulling the strap on the underside and closing the buckle.



EN

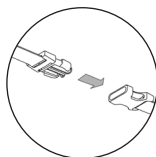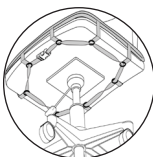On a four-legged chair or stool, the cushion is attached by using the cord locks or by tying a knot.

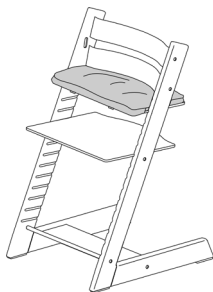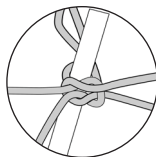
On a high chair, the cushion is attached by closing the two buckles on the underside of the plate.

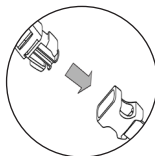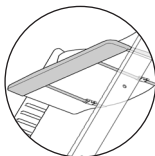
Full foot support is important when sitting on the cushion. In case the child can't reach the ground, use a stable foot support.

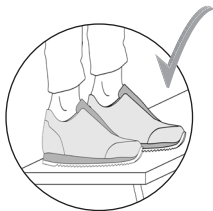During normal use of the cushion, the child will move around in the seat, and sitting on the balls will therefore not feel uncomfortable.

The frame of the chair is not an optimum foot support.

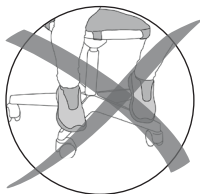Use environment

Protac Ball Cushion® is for indoor use. The cushion can be used both in institutions and private homes.

Indications

Protac Ball Cushion® is a sensory-stimulating aid intended to alleviate difficulties concentrating and help a seated user train balance. Furthermore, it can be used for persons (children, adults, elderly) with sensory impairment.

Protac Ball Cushion® contains balls which move and apply dynamic pressure on to the body, thereby stimulating the muscle-joint sense and the tactile sense. This increases body awareness and calms the brain and body, and the product thereby alleviates physical and psychological unrest.

Contra-indications

Persons with pressure ulcers or persons at risk of developing pressure ulcers are not to use the product.

Not suitable for wheelchair users.

Not for unsupervised use by users with paralysis and poor sensitivity.

Maintenance

It is recommended to air the cushion on a regular basis.

Cleaning

Protac Ball Cushion® can be washed in a regular household washing machine. The frequency of washing depends on use. When washing, please follow the instructions on the care label. Washing instructions can also be found on page 51 of this manual.

Protac does not recommend the use of stain remover.

Key to symbols - care label

EN

	Maximum washing temperature		Do not wash
	Can be tumble dried		Do not tumble dry
	Dry flat		Do not iron
	Line dry in the shade		Do not dryclean
	Shrinkage to be expected		Do not bleach
	Minimum capacity 7 kg		Minimum capacity 16 kg

Service life and warranty

The expected service life of the product is 5 years. This is based on correct use, cleaning, washing and maintenance.

Products used beyond the expected service life must regularly undergo a maintenance inspection which implies checking for wear and tear in fabric and seams.

The product comes with a 2-year warranty, see page 2.

The warranty will be invalid in case of:

- overload or excessive use of the product
- lack of or incorrect maintenance of the product

Disposal

The product can be broken down into the following components for disposal:

- The plastic balls can be disposed of as hard plastic
- The fabric can be disposed of as textiles

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Protac A/S